



Rec'd PCT/PTO 01 APR 2005
PCT/EP 03 / 11360

10/530073

INVESTOR IN PEOPLE

PRIORITY DOCUMENT
SUBMITTED OR TRANSMITTED IN
COMPLIANCE WITH
RULE 17.1(a) OR (b)

REC'D 06 NOV 2003

WIPO


The Patent Office
Concept House
PCT
Cardiff Road
Newport
South Wales
NP10 8QQ

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

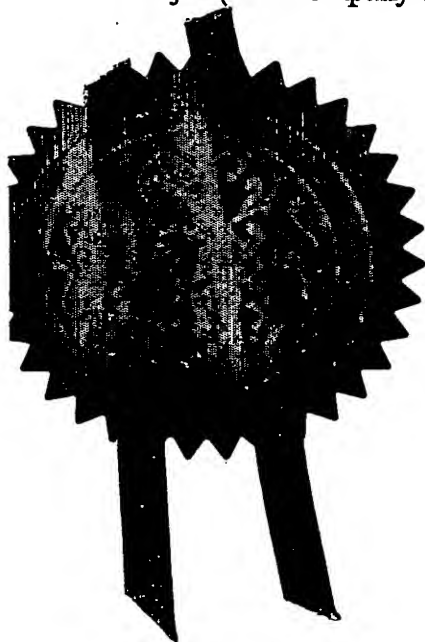


Signed *A. B. Jones*

Dated 1 October 2003

BEST AVAILABLE COPY

Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

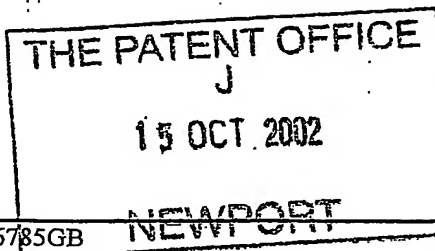

The Patent Office

Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

1. Your reference	WYET / P26785GB		
2. Patent application number (The Patent Office will fill in this part)	15 OCT 2002	0223960.6	16OCT02 E755930-1 D02866 P01/7700 0.00-0223960.6
3. Full name, address and postcode of the or of each applicant (underline all surnames)	Wyeth Five Giralda Farms Madison New Jersey 07940 0874 United States of America		
Patents ADP number (if you know it)	08363970001		
If the applicant is a corporate body, give the country/state of its incorporation	United States of America		
4. Title of the invention	ASSAY METHODS		
5. Name of your agent (if you have one)	ERIC POTTER CLARKSON		
"Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)	PARK VIEW HOUSE 58 THE ROPEWALK NOTTINGHAM NG1 5DD		
Patents ADP number (if you know it)	1305010		
6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (if you know it) the or each application number	Country	Priority application number (if you know it)	Date of filing (day / month / year)
7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application	Number of earlier application		Date of filing (day / month / year)
8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (Answer 'Yes' if:	YES		
a) any applicant named in part 3 is not an inventor; or b) there is an inventor who is not named as an applicant, or c) any named applicant is a corporate body. See note (d))			

Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form.
Do not count copies of the same document

Continuation sheets of this form

Description 29

Claims(s) 4

Abstract 1

Drawing(s) 13

10. If you are also filing in any of the following, state how many against each item.

Priority Documents 0

Translations of priority documents 0

Statement of inventorship and right to grant of a patent (Patents Form 7/77) NO

Request for preliminary examination and search (Patents Form 9/77) NO

Request for substantive examination (Patents Form 10/77) NO

Any other documents (please specify)

11. I/We request the grant of a patent on the basis of this application.

Eric Potter Clarkson

Signature
ERIC POTTER CLARKSON

Date
14 October 2002

12. Name and daytime telephone number of person to contact in the United Kingdom 0115 9552211

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 01645 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

ASSAY METHODS

The present invention relates to assay methods and, in particular, to methods for detecting the presence of a virus, especially Marek's disease virus (MDV), in an avian tissue sample.

Background

Marek's disease virus (MDV) is a herpesvirus, which causes lymphoproliferative disease in chickens. Even after the introduction of vaccines against MDV, the infection still causes considerable losses in the poultry industry. MDV is divided into three serotypes, all of which establish latent infections. Serotype 1 includes oncogenic viruses, serotype 2 non-oncogenic viruses and serotype 3 includes the turkey herpesviruses (HVT) (Bülow *et al* (1976) *Zentralblatt für Veterinarmedizin*, 23B, 391-402).

The traditional diagnosis of Marek's disease is based on the clinical signs and pathological alterations. However, more specific methods for surveillance of the prevalence of MDV would be desirable. The detection of viral antigen in the feather follicle epithelium by the agar gel precipitation test has been described by Haider *et al* (1970) *Poultry Science*, 49, 1654-1657. The different serotypes can be differentiated by the agar gel precipitation test (Lee *et al* (1983) *Journal of Immunology*, 130, 1003-1006), but the sensitivity of that test is inferior to that of enzyme-linked immunosorbent assay (ELISA) and DNA hybridization (Davidson *et al* (1986) In: Current research on Marek's disease. *Proceedings of the 5th International Symposium on Marek's Disease* (pp. 311-316). Tallahassee: Rose Printing Company, Inc.).

The preferred samples for virus isolation are buffy-coat cells, which can be co-cultivated with susceptible primary cell cultures. Immunofluorescent assay (Kitamoto *et al* (1979). *Biken Journal*, 4, 137-142) or ELISA (Cheng *et al* (1984) *Avian Diseases*, 4, 900-911), can be used for subsequent identification of the MDV serotype. Alternatively, the serotype can be identified by restriction endonuclease analysis (Ross *et al* (1983). *Journal of General Virology*, 64, 2785-2790) or polymerase chain reaction (PCR) (Wang *et al* (1993) *Molecular and Cellular Probes*, 7, 127-131. *In situ* hybridization has been used for detection of MDV genome in infected tissue (Endoh *et al* (1996) *Journal of Veterinary Medical Science*, 58, 969-976; Ross *et al* (1997) *Journal of General Virology*, 78, 2191-2198), but this technique is probably too laborious for routine diagnoses. (Davidson *et al* (1995) *Avian Pathology*, 24, 69-94; and Davidson *et al* (1996) In: Current research on Marek's disease. *Proceedings of the 5th International Symposium on Marek's Disease* (pp. 311-316). Tallahassee: Rose Printing Company, Inc.), applied MDV serotype 1-specific PCR techniques to full blood and tumour tissue samples from commercial chicken and turkey flocks, the majority of which had neoplastic disease. Wang *et al* (1993) *Molecular and Cellular Probes*, 7, 127-131; Young, P. & Gravel, J. (1996) In Current research on Marek's disease. *Proceedings of the 5th International Symposium on Marek's Disease* (pp. 308-310),. Tallahassee: Rose Printing Company, Inc.; and Silva, R.F. & Witter R.L. (1996) In Current research on Marek's disease. *Proceedings of the 5th International Symposium on Marek's Disease* (pp. 302-307). Tallahassee: Rose Printing Company, Inc., applied a MDV serotype 1-specific PCR protocol to various tissues of chickens experimentally inoculated with the JM/102 strain:

Handberg *et al* (2001) *Avian Pathology* 30 : 243-249 describe the use of serotype 1- and serotype 3-specific PCR for the detection of MDV in

chickens. Tissue samples were taken from blood (buffy-coat cells), spleen, liver, skin, feather tips and ovaries.

Description of the invention

5

The present invention provides further methods for detecting a virus, especially MDV, in avian tissue samples.

10 In a first aspect the invention provides a method of detecting a virus in an avian tissue sample comprising: extracting genetic material from an avian tissue sample; and testing the extracted genetic material to detect any genetic material from the virus; characterised in that the avian tissue sample is derived from one or more feathers of the axillary tract.

15 By using feather samples the test can be carried out on live animals. Sampling of feathers is simple, quick and practical under field conditions. Feather samples can be placed in a suitable container and tested immediately or stored for future testing, as desired. By contrast, sampling blood requires great care to prevent blood clots forming, including transport
20 of blood under cool, controlled conditions. Blood clotting leads to negative test results. Internal organ samples, such as spleen and tumour samples, must be transported on wet ice, which is impractical under field conditions.

25 By selecting axillary tract feathers from which to derive a tissue sample, the invention provides significant advantages over known methods which take tissue samples from different parts of the bird.

30 Surprisingly, virus can be detected in axillary tract feathers according to the invention when it cannot be detected in other tissue samples, including other feathers. Accordingly, the methods of the invention are particularly suitable

for monitoring the extent to which a flock of birds has been immunised effectively with MDV vaccine, by detecting the presence of the vaccine strain in axillary tract feather tissue samples.

- 5 By "avian" we include any bird, but preferably birds which are produced commercially, especially poultry such as chickens, turkeys, ducks, etc.

By "axillary tract feathers", we include the meaning of the feathers located in the region of a bird marked "axillary" in the accompanying feathering diagram (Figure 1). Preferably, the axillary tract feather selected is a "pin feather", that is, an immature growing feather. The term "pin feather" will be familiar to skilled persons. For example, van Tyne J & Berger AJ (1959) *Fundamentals of Ornithology*, John Wiley, New York refer to a pin feather as "a new, growing feather, still not completely unsheathed". Lucas AM & Stettenheim PR (1972) *Avian Anatomy*, Integument Part 1, US Government Printing Office, Washington, pp 199-200 remark "The new feather is tightly furled inside a sheath while it forms. As it appears above the skin, it has a long conical shape with a blunt tip and a slightly moist surface. A feather at this stage in any generation is often called a pin feather". (See also Figure 1(a), which is taken from Lucas & Stettenheim (1972)).

Preferably, the axillary tract feathers are taken from chicks which are advantageously less than a month old.

- 25 Preferably, the method provides quantitative information on the amount of virus, especially MDV, in the sample.

Preferably, the method is specific for MDV serotype 1, and more preferably the method is specific for MDV-1 Rispens strain CVI 988. The latter strain is a commercial vaccine strain produced by Fort Dodge, IA, USA which is

available from the American Type Culture Collection (ATCC), Mannassas, VA, USA.

Advantageously, the method involves the use of a PCR reaction.

5 Preferably, before said PCR reaction is carried out, the extracted genetic material to be tested is treated with an agent to overcome the inhibitory effect of any feather tissue factor which may be present. This inhibitory effect appears to be associated with melanin and can therefore be a particular problem when feathers from brown birds are sampled. Preferably,
10 the agent is selected from one or more of bovine serum albumin; porcine (pig) albumin; and ovine (sheep) albumin.

Skilled persons will be aware of a range of detection methods for detecting viral, especially MDV, genetic material which could be used in the methods
15 of the invention, such as the methods of Handberg *et al* (2001) *supra*. A particularly preferred method for detecting MDV-1 strain CVI 988 is as follows:

- 20 (i) providing forward and reverse primers for a nucleic acid polymerase, which primers are selected from the nucleotide sequence which flanks the 132 bp repeat nucleotide sequence of MDV;
- (ii) amplifying nucleic acid sequences between the primers;
- 25 (iii) detecting the number of 132 bp repeat sequences in the amplified nucleic acid sequences; and
- (iv) relating the number of 132 bp repeat sequences to the identity of the viral nucleic acid and thereby identifying the type of MDV in the tissue sample, multiple copies of the 132 bp repeat sequence being indicative of MDV-1 strain CVI 988.

A preferred quantitative method for use in detecting MDV according to the present invention comprises:

- 5 (a) providing a polynucleotide sequence which is capable of binding specifically to a MDV-specific target polynucleotide;
 - (b) contacting the extracted genetic material with a probe whereby the probe binds specifically to its target MDV polynucleotide;
 - 10 (c) determining whether the probe has bound to its target MDV polynucleotide; and
 - (d) determining whether the sample contains MDV on the basis that the presence of the target polynucleotide indicates the presence of MDV in the sample.
- 15 Step (d) preferably provides a quantitative determination of the amount of virus in the sample.

Advantageously the step of determining whether the probe has bound to a target polynucleotide comprises amplifying a region of the target
20 polynucleotide, which region comprises the binding site of the probe.

Preferably the probe has the sequence 5' AGA CCC TGA TGA TCC GCA TTG CGA CT 3'.

25 Preferably, amplification is primed by the following primers:

Forward primer (GGT CTG GTG GTT TCC AGG TGA) which is located at NT positions 1341-1361 in the GA strain Meq gene sequence. The GA (Georgia) strain was a 1964 isolate from Georgia, from an
30 ovarian tumour. Reference: C.S. Eidson & S.C. Schmittle (1968).

Studies on acute Marek's disease. I. Characteristics of isolate GA in chickens. Avian Diseases 12;467-476.

Reverse primer (GCA TAG ACG ATG TGC TGC TGA) is located at
5 NT positions 1413-1393.

Advantageously, the probe is labelled fluorescently and the step of determining whether the probe has bound to a target polynucleotide comprises determining the fluorescent emissions of the probe.

10

Two fluorescent dyes are brought into physical proximity by direct conjugation at opposite ends of a short oligo probe (5' reporter fluorochrome, usually 6-FAM, and 3' quencher fluorochrome, usually TAMRA). When the high-energy fluorophore (FAM) is excited at 488 nm, instead of the expected fluorescence emission at 520 nm the captured
15 energy is transferred to the lower energy fluorophore (TAMRA) and is emitted at 580 nm (fluorescence resonance energy transfer, FRET, has occurred). Using the fluorescein/rhodamine reporter/quencher combination, FRET will effectively occur even when the groups are separated by 25-30
20 bases of DNA. During the course of a TaqManTM assay the two fluorophores are physically detached from each other by the 5'-nuclease action of Taq DNA polymerase - after which 488 nm stimulation results in visible FAM emission at 520 nm.

25 Dual-labeled probes usually have a 5'-reporter dye, such as FAM, TET, HEX, JOE or VIC and a 3'-quencher group, such as TAMRA or Dabcyl (a universal quencher).

Suitable fluorescent labels are within the common general knowledge of
30 skilled persons. The following reagents are available from Sigma-Aldrich, UK: HEX stands for hexachloro-fluorescein; TET stands for tetrachloro-

fluorescein; Joe is 6-carboxy-4',5'-dichloro-2',7'-dimethoxyfluorescein; ROX is 5-carboxy-Rhodamine; dabcyl is 4-((4-(dimethylamino)phenyl)azo)benzoic acid.

5 Structures and further information on all of these are available on the molecular probe website www.probes.com. VIC dye is available from Applied Biosystems.

6-FAM = 6-carboxyfluorescein; (a phosphoramidite); yellow-green dye;
10 absorbance maximum = 494 nm, emission maximum = 525nm

HEX = a phosphoramidite; pink dye; absorbance maximum = 535, emission maximum = 556

TET = a phosphoramidite; orange dye; absorbance maximum = 521,
15 emission maximum = 536

VIC = absorbance maximum 538, emission maximum = 554

JOE = absorbance maximum 521, emission maximum = 547

TAMRA = 6-carboxy-tetramethyl-rhodamine; absorbance maximum = 555,
emission maximum = 580

20 Dabcyl = absorbance maximum = 453, no maximum emission (universal quencher)

From the foregoing description it will be apparent that an important aspect of the methods of the invention is the use of avian tissue samples which are
25 more convenient and useful than avian tissue samples used previously for detecting viral, especially MDV, infection. Accordingly, further aspects of the invention relate to the provision of avian tissue samples.

In a second aspect the invention provides an isolated avian tissue sample
30 from one or more feathers from the axillary tract.

By "isolated" we include the meaning that the tissue sample is free of a substantial amount of the material with which it is normally associated in nature. For example, the tissue sample may be stored in a container, or
5 be derived from the original axillary tract feather by a variety of isolation; extraction and/or purification methods.

Preferably the isolated tissue sample consists of the proximal portion (the non-barbed portion which is attached to the skin and which contains the
10 pulp -- see accompanying figures). Accordingly it is preferred that the proximal portion of the axillary tract feather is isolated from the distal (barbed) portion of the feather. This is easily achieved simply by cutting off the proximal portion of the feather with a pair of scissors and discarding the distal portion.

15

In a third aspect the invention provides a genetic material-containing extract from an avian tissue sample wherein the extract is taken from a sample of tissue as described in relation to the second aspect of the invention.

20 It will be appreciated that samples according to the second and third aspects of the invention may be collected and/or prepared in the field, or transported to a separate location, such as a laboratory, for preparation and testing. Hence, further aspects of the invention relate to samples according to the second and/or third aspects of the invention stored in a form suitable for
25 transport to a separate location.

The feathers could be stored complete (eg. in 20 ml Sterilin universal plastic tubes), or after cutting off the proximal portion required for DNA preparation (eg. in 1 ml Eppendorf snap-cap or screw-cap tubes).
30 Alternatively, the feathers could be stored in heat-sealed or tied plastic bags.

For short-term storage (eg. 1 week), the feathers could be stored at 4°C, but, for longer periods of storage, they should be stored frozen at -20°C.

The DNA should be stored at -20°C in screw-cap or snap-cap 0.5 ml or 1.5 ml Eppendorf tubes, in Tris-EDTA buffer or, if the DNA is to be used in
5 Taqman analysis, in water since EDTA will inhibit the Taqman reaction.

Advantageously, the results of the methods are furnished in an intelligible format. Preferably, the results are recorded or stored on an information carrier. However, the step of furnishing the results could be by
10 communicating the results orally.

By "information carrier", we include any means of storing information, such as paper, a computer disk; an internet-based information transfer system, such as an e-mail or internet page, or electronic file, etc. Of course,
15 an "intelligible format" is also intended to embrace encrypted information which can be deciphered with an approximate key.

Specific description and drawings/figures

20 Examples embodying certain aspects of the invention will now be described with reference to the following figures in which:

Figure 1(a) shows the development of feathers above the skin, during the first four generations. A circle indicates a feather that is still growing; a dot,
25 a feather that is fully grown.

Figure 1(b) is a diagram of the feathering pattern in chickens which shows the axillary tract.

The axillary tracts lie on the underside of the chick and extend, each side, from the lower neck to the upper abdomen, underneath the wings and alongside the breastbone. Figure 1 is taken from the book 'Bird Structure: An approach through evolution development and function in the fowl'.
5 D.A.Ede, Publisher: Hutchinson Educational, 1964. The feathering pattern is the same in all birds.

Figures 2(a) to 2(d) are photographs of the axillary tract feathers and close ups of individual feathers showing the portions which are retained to
10 provide tissue samples according to the invention. Figures 2(a) and 2(b) demonstrate the axillary tract feathers. Figures 2(c) and 2(d) demonstrate the part of the feather taken for analysis'.

Figure 3:

15 This experiment was performed to test/optimize the use of the *Meq* Taqman primer/probe set. BAC10 DNA was used as the target DNA since it was known to contain many copies of the viral genome.

Taqman assay was used to detect the *Meq* gene in DNA derived from Rispens virus genome cloned into a Bacterial Artificial Chromosome
20 (BAC10). The DNA was used in ten-fold dilutions. During each cycle of real-time PCR, the reporter fluorochrome FAM is released and able to fluoresce. Therefore, with each cycle, fluorescence intensity increases. The Ct value (the cycle at which fluorescence passes a fixed threshold) is a measure of the starting copy number of the target sequence. The lower the
25 Ct value, the higher the starting copy number of target sequence. A Ct value of 40 indicates that no target sequence, or an undetectable amount of target sequence, is present. This figure shows that dilutions of this DNA preparation between 1:10 and 1:10000 gave a detectable amount of *Meq* PCR product. The Ct value increases linearly with increasingly dilute DNA.

30

Figure 4:

This experiment was performed to test/optimize the use of the *Meq* Taqman primer/probe set on DNA derived from Rispens- MDV infected cells.

Taqman assay was used to detect the *Meq* gene in DNA derived from BAC10, and in DNA derived from Rispens MDV-infected chick embryo
5 fibroblast (CEF) cells. The DNA was used in ten-fold dilutions. This figure shows that for both BAC10 DNA and Rispens-infected CEF DNA, the Ct value increases linearly with increasingly dilute DNA.

Figure 5:

10 This experiment was performed to test/optimize the use of the *Meq* Taqman primer/probe set on DNA derived from tissue samples from Rispens- MDV infected chickens.

Taqman assay was used to detect the *Meq* gene in DNA from a spleen of a Rispens-inoculated chick (11 dpi), and an age-matched uninoculated chick.
15 DNA from Rispens-infected CEF cells was used as a positive control. The DNA was used in ten-fold dilutions of 1 mg/ml stocks and 1 µl was used per reaction. Although the Ct values for the uninoculated spleen DNA were lower than 40, they clearly did not increase with increasing concentrations of DNA. However, *Meq* detection in the inoculated spleen DNA rose
20 significantly above this baseline when the DNA was used neat, or at 1:10 dilution. We thus established that, for Taqman analysis of DNA taken from tissue samples of MDV-inoculated chicks, we would use 1 µl DNA from a 1 mg/ml stock (ie. 1 µg) DNA per reaction.

25 *Figure 6:*

This experiment was performed to follow the time-course of Rispens MDV-infection in feather axillary tracts of inoculated chicks by Taqman assay.

Taqman assay was used to detect the *Meq* gene in DNA prepared from
30 feather tips of chicks at 0 (uninoculated), 10, 15, 20 and 28 dpi post inoculation with Rispens. A group of five chicks were sampled at each

time-point (four chicks at 0 dpi). 1 µg DNA was used in the Taqman assay. Mean Ct values for each group are plotted. The Ct values decrease from 0 – 15 dpi, then increase again from 20 – 28 dpi, showing a peak of infection 15 – 20 dpi.

5

Figure 7:

This experiment was performed to follow the time-course of Rispens MDV-
10 infection in feather axillary tracts of inoculated chicks by 132 bp repeat PCR.

DNA was prepared from feather tips of Rispens-inoculated chicks at 0 (uninoculated), 10, 15, 20 and 28 dpi. A group of five chicks were sampled at each time-point (four chicks at 0 dpi). 1 µg DNA was used in PCR
15 (including 10 µg BSA per reaction) and the samples run on a 1% agarose gel containing ethidium bromide.

(M) = Lambda molecular size markers, (-) = water (negative control), (+) = Rispens BAC10 DNA (positive control). Days post infection are indicated underneath the gels.

20 (a) To confirm the PCR-quality and quantity of each DNA sample, PCR was performed to detect an endogenous retrovirus sequence present in all chicken cells. The 360 bp endogenous retrovirus product was detected in all of the feather samples confirming the PCR-quality of each samples.

(b) 132 bp repeat PCR was performed. The lane marked (+) shows the 132
25 bp repeat ladder PCR products obtained with Rispens BAC10 – six copies are clearly distinguished. The negative control shows no PCR product, as do three of the four uninoculated chicks. The fourth uninoculated chick shows a faint product band equivalent to 3 copies of the repeat, indicating that the chick was contact-infected by the inoculated chicks housed in the
30 same room. All of the inoculated chicks were 132 bp repeat positive at 11, 15, 20 and 28 dpi. The PCR product representing a certain repeat copy

number predominated in many cases and the number of copies represented by this predominant product varied between samples. This indicates that sub-clones of the inoculum virus, with a set number of repeats, come to predominate in different chicks.

5

Figure 8:

Gives the nucleotide and amino acid sequences of the *meq* gene of MDV serotype 1. This sequence is described by Jones *et al* (May 1992) *Proc. Natl. Acad. Sci USA*, vol 89, pp 4042-4046 and has been deposited with
10 GenBank under accession no. M89471.

EXAMPLE 1: quantitative PCR assay for MDV-1

15 TaqmanTM quantitative PCR is an established technique used to quantify the amount of starting PCR target by determining the number of PCR cycles required to reach a fluorescence threshold (defined mathematically by the Ct value). A higher copy number of target sequence in the sample requires fewer PCR cycles to reach the Ct threshold.

20

The Taqman primers and probe were designed from the *Meq* gene sequence of the MDV strain GA (see later). This sequence is published in: Jones D, Lee L, Liu JL, Kung HJ, Tillotson JK. Marek disease virus encodes a basic-leucine zipper gene resembling the fos/jun oncogenes that is highly
25 expressed in lymphoblastoid tumors. *Proc Natl Acad Sci U S A*. 1992 May 1;89(9):4042-6. The Sequence Accession No. is: M89471. The Applied Biosystems 'Primer Express' software was used to select the optimum primer/probe sequences from the *Meq* sequence.

The specific primers used in this analysis multiply a 73bp sequence in the meq gene of MDV serotype 1 virus. This region is common to both vaccine strains and field isolates.

- 5 The experiments described use CVI 988 vaccine as an example of serotype 1 Marek's disease virus. The results can be applied to all serotype 1 Marek's viruses because the area of the sequence the primers are directed against is conserved in MDV serotype 1 viruses.

10 *Protocol for TaqmanTM Meq gene PCR analysis on feather tip samples*

1. Materials required:

Reagents:

- 15 • TNE buffer (store at room temperature) contains Tris (10 mM), NaCl (150 mM), EDTA (1 mM) and the pH is adjusted to pH 7.5 using HCl
- Sodium Dodecyl Sulphate (SDS) 10% solution (store at room temperature)
- Proteinase K (lyophilised powder from Sigma # P-6556, and make up a stock of 20 mg/ml in water, stored at -20°C)
- 20 • Phenol pH 7.9 obtained from Sigma (catalogue no. P-4557) stored at 4°C)
- Chloroform (stored at -20°C)
- 3M Sodium Acetate pH5.2 (stored at room temperature)
- 25 • Filtered neat ethanol (stored at room temperature)
- Filtered 70% cold ethanol (stored at 4°C)
- PCR quality water
- PCR quality water containing 800 µg/ml Bovine Serum Albumin (BSA), filtered
- 30 • Ice

- TaqmanTM PCR core reagents – TaqmanTM buffer, MgCl₂, dNTPs, Taq polymerase, Uracil N-Glycosylase (Perkin Elmer Biosystems)
Meq forward primer 5' GGT CTG GTG GTT TCC AGG TGA 3'
(MWG Biotech)
- 5 Meq reverse primer 5' GCA TAG ACG ATG TGC TGC TGA 3'
(MWG Biotech)
Meq probe 5' FAM AGA CCC TGA TGA TCC GCA TTG CGA CT 3'
TAMRA
(Sigma-Genosys Ltd)
- 10 (FAM & TAMRA are the fluorescent tags)

- 6-FAM = 6-carboxyfluorescein; (a phosphoramidite); yellow-green dye;
absorbance maximum = 494 nm, emission maximum = 525nm
TAMRA = 6-carboxy-tetramethyl-rhodamine; absorbance maximum = 555,
15 emission maximum = 580

- 5'FAM-3'TAMRA labelled probes are available from: Sigma-Genosys Ltd.
(London Road, Pampisford, Cambridgeshire, CB2 4EF, UK. Tel. 01223
839200)
- 20 5'VIC-3'TAMRA labelled probes are available from: Applied Biosystems
Ltd. (Kelvin Close, Birchwood Science Park North, Warrington, Cheshire,
WA3 7PB)

- 25 Other suppliers are:
 - Integrated DNA Technologies (IDT) 1710 Commercial Park, Coralville,
IA, 52241, USA
 - Oswel Research Products Ltd.: Lab 5005, Medical and Biological
Sciences Building, University of Southampton, Boldrewood, Bassett
30 Crescent East, Southampton, SO16 7PX, Tel: 02380 592984

Equipment

- Sterilin 20 ml plastic universal tubes
- Clean scissors & forceps
- 5 • Water bath set to 50°C
- Micro centrifuge
- 1.5 ml snap cap Eppendorf tubes (autoclaved)
- 1.5 ml screw-cap Eppendorf tubes (autoclaved)
- 0.5 ml snap-cap Eppendorf tubes (autoclaved)
- 10 • Vacuum/freeze-drier (not essential)
- Spectrophotometer
- Dedicated PCR cabinet, pipettes and autoclaved tips
- Thermo-fast 96-well PCR plate and caps (Perkin Elmer/Applied Biosystems)
- 15 • ABI Prism 7700 Sequence Detector (Perkin Elmer/Applied Biosystems)

2. Collecting the feathers:

- Pluck 8 – 10 'pin' feathers (short, newly growing feathers with plenty of pulp) from the brachial feather tract of each chicken (see figures).
- 20 • Place feathers in a plastic 'universal' tube for transport back to the laboratory.

3. DNA preparation from feather tips:

- Cut off and save the proximal 1 cm of the feather (i.e. the non-barbed part which is attached to the skin and which contains the pulp – see photographs). Discard the distal barbed part of the feather.
- 25 • For each chicken, place the 8 – 10 saved feather ends in a 1.5 ml snap-cap Eppendorf tube.

Add 500 µl proteinase K sample buffer (TNE buffer containing 0.5% SDS) containing 100 µg proteinase K (add proteinase K just before use).

30

- Incubate at 50°C in a water-bath for 1.5 – 2 hours.
- Microcentrifuge the tubes (6000 rpm, 10 min), to 'pellet' feather tips & debris.
- Transfer supernatant to a new snap-cap tube (if the feathers were from brown birds, the supernatants will be brown due to the presence of melanin).
- Add an equal volume (500 µl) of phenol to the supernatant again.
- Vortex
- Centrifuge at 13000 rpm, 2 min.
- Transfer the upper (aqueous) phase to a new snap-cap tube.
- Add an equal volume (500 µl) of phenol to the supernatant again.
- Vortex
- Centrifuge at 13000 rpm, 2 min.
- Transfer the upper (aqueous) phase to a new snap-cap tube.
- Add an equal volume (500 µl) of cold chloroform.
- Vortex.
- Centrifuge at 13000 rpm, 2 min.
- Transfer the upper phase to a 1.5 ml screw-cap tube.
- Add 1 ml filtered 100% ethanol.
- Add 50 µl of 3M Sodium Acetate.
- Gently mix by inverting the tube, and leave at room temperature for 20 minutes (the DNA will become visible as it precipitates).
- Centrifuge 13000 rpm, 2 minutes, to pellet the DNA (if white chickens were used, pellet will be white; if brown chickens used, pellet brown).
- Discard the supernatant.
- Rinse the pellet twice with 500 µl of 70% cold ethanol, by gently running the ethanol down the side of the tube, then pouring off (take care not to dislodge the DNA pellet).

- Cover the open top of the tube with Parafilm, and make several puncture holes in the Parafilm using a needle.
- Place tubes in a vacuum drier for about 5 minutes to dry the pellet (alternatively air-dry).

5

Re-suspend the pellet in 50 μ l PCR quality water by gently vortexing.

Determine the concentration of the DNA preparation using a spectrophotometer.

Adjust the concentration to 1 mg/ml in water.

10 Store at -20°C .

4. *TaqMan*TM quantitative PCR assay (Perkin Elmer Biosystems)

- Set up reactions in a PCR-dedicated cabinet, using PCR-dedicated pipettes and autoclaved tips
- 15 • Work on ice
- Prepare master mix containing the following reagents for the appropriate number of samples – set up duplicate reactions for each sample (volumes given per reaction):

Component	Volume per 25 μ l reaction	Final Concentration
Taqman buffer	2.5 μ l	
MgCl ₂ (5 mM)	5.0 μ l	1.0 mM
dATP (10mM)	0.5 μ l	0.2 mM
dCTP (10 mM)	0.5 μ l	0.2 mM
dGTP (10mM)	0.5 μ l	0.2 mM
dUTP (10mM)	0.5 μ l	0.2 mM
Water containing 800 μ g/ml BSA	11.6 μ l	~10 μ g BSA/reaction

Meq probe (10 μ M)	0.5 μ l	0.2 μ M
Meq forward primer (10 μ M)	1.0 μ l	0.4 μ M
Meq reverse primer (10 μ M)	1.0 μ l	0.4 μ M
Taq Gold DNA pol (5U/ μ l)	0.13 μ l	26 U/ml
Uracil N-glycosylase (1U/ μ l)	0.25	10 U/ml

- Vortex to ensure complete mixing
- Place a thermo-fast PCR plate into a plate holder on ice and add 24 μ l master mix to each well to be used
- Add 1 μ l autoclaved water to no-template-control (NTC) wells and cap these wells prior to opening any DNA samples
- Add 1 μ l positive control DNA (= 1 μ l of 1 mg/ml preparation), eg. DNA from Rispens-infected CEF to appropriate wells and cap these wells
- Add 1 μ g sample DNA (= 1 μ l of 1 mg/ml preparation), to appropriate wells and cap
- Briefly pulse plate in centrifuge
- Place plate in ABI Prism 7700 Sequence Detector (Applied Biosystems) and set up computer to read FAM fluorescence (Meq probe), run samples
- *Thermocycling parameters*
 - 50°C 2 min
 - 95°C 10 min
 - 94°C 15 sec) x 40
 - 60°C 1 min)

- Analyse data using Microsoft excel – for each sample there will be a Ct value (the PCR cycle at which the amount of fluorescent product is first detected above baseline level); calculate mean Ct value for duplicates
5 for each DNA sample.

Results of the above experiment are shown in Figure 6.

10 Notes

The Forward primer (GGT CTG GTG GTT TCC AGG TGA) is located at NT positions 1341-1361 in the GA strain Meq gene sequence. The reverse

15 primer (GCA TAG ACG ATG TGC TGC TGA) is located at NT positions 1413-1393. The probe was designed to specifically anneal between the two primers on the Meq target sequence (see Figure 8).

During PCR, the fluorogenic probe binds between the two primers and, during each extension cycle, the 5' nuclease activity of the Taq polymerase
20 cleaves the probe, releasing the reporter fluorochrome FAM that is then able to fluoresce. Therefore, with each cycle, fluorescence intensity increases. The Ct value (the cycle at which fluorescence passes a fixed threshold) is a measure of the starting copy number of the target sequence: the higher the starting copy number, the lower the Ct value.

25

Use of Bovine Serum Albumin during PCR: feather tissues (especially those from brown chickens) contain melanin, which has been shown to be inhibitory to PCR. Use of BSA in the reaction overcomes this melanin-induced inhibition. Experiments were carried out in brown chickens. The

addition of BSA followed the method of Giambernardi *et al* (1998)
Biotechniques 25: 564-6.

**EXAMPLE 2: Method for detecting specific MDV-1 (CV1 988 vaccine)
strain**

5 This example of the invention relates to the use of a specially modified version of the 132 base pair (bp) repeat polymerase chain reaction (PCR) test to detect the presence of CVI 988 Marek's vaccine in chicken feathers.

10 The 132 base pair repeat genetic sequence is located in the internal repeat long (IR₁) segment of the Marek's disease virus (serotype 1) genome. The complete genomic sequence of MDV 1 is described in Tulman *et al.* (Sept 2000) *J Virol.* Vol. 74 No. 17, p7980-7988 and has been deposited in GenBank under accession no. AF 243438. CV1 988 isolates (vaccine strains) of serotype 1 Marek's disease have multiple copies of this repeat
15 segment, whilst field strains have single copies (Silva *et al* (1992) *Avian Dis* 36 : 521-528; and Becker *et al* (1993) *Virus Genes* 7 : 277-287). Measuring the number of copies affords the possibility of differentiating vaccine strains from field strains. (Becker *et al* (1992) *J Virol Methods* 40 : 307-322 and Kopacek *et al* (1993) *Acta Virol* 37 : 191-195).

20

Feathers are sampled according to the enclosed figures, the proximal tips of the axillary tract feathers being used. PCR analysis (see below) demonstrated multiple copies of the 132 bp segment in animals between 11 and 28 days post vaccination. Vaccination was carried out on birds at one
25 day of age.

An acceptable variation of the test is to use feathers in birds of any age either taken as fresh samples, or stored for testing at a later date to determine the presence of the CV1988 vaccine.

Protocol for 132 bp repeat PCR on Feather tip samples

1, Materials required:

5

Reagents:

- TNE buffer (store at room temperature) contains Tris (10 mM), NaCl (150 mM), EDTA (1 mM) and the pH is adjusted to pH 7.5 using HCl
- Sodium Dodecyl Sulphate (SDS) 10% solution (store at room
10 temperature)
- Proteinase K (lyophilised powder from Sigma # P-6556, and make up a stock of 20 mg/ml in water, stored at -20°C)
- Phenol pH 7.9 obtained from Sigma (catalogue no. P-4557) (stored at 4°C)
- 15 • Chloroform (stored at -20°C)
- 3M Sodium Acetate pH5.2 (stored at room temperature)
- Filtered neat ethanol (stored at room temperature)
- Filtered 70% cold ethanol (stored at 4°C)
- PCR quality water
- 20 • PCR quality water containing 800 µg/ml Bovine Serum Albumin (BSA), filtered
- Taq gold DNA polymerase (5 U/µl), Taq buffer, MgCl₂ (25 mM) from Bio/Gene Ltd, Kimbolton, Cambridgeshire, England
- dATP, dTTP, dCTP, dGTP (100 mM stocks) obtained from Promega
25 (USA); we prepare a mix containing all four of these at 10 mM each, stored at -20°C
- DNA molecular size markers
- Agarose and TBE buffer

Primer sequences:

MD-132 FOR; 5' TACTTCCTATATAGATTGAGACGT-3'

MD-132 REV: 5' GAGATCCTCGTAAGGTGTAATATA-3'

5

Equipment:

- Sterilin 20 ml plastic universal tubes
- Clean scissors & forceps
- Water bath set to 50°C
- 10 • Micro centrifuge
- 1.5 ml snap cap Eppendorf tubes (autoclaved)
- 1.5 ml screw-cap Eppendorf tubes (autoclaved)
- 0.5 ml snap-cap Eppendorf tubes (autoclaved)
- Vacuum/freeze-drier (not essential)
- 15 • Spectrophotometer
- Dedicated PCR cabinet, pipettes & autoclaved tips
- Thermal cycler
- Agarose gel apparatus

20

2. Collecting the feathers:

- Pluck 8 – 10 'pin' feathers (short, newly growing feathers with plenty of pulp) from the axillary feather tract of each chicken (see figures).
- Place feathers in a plastic 'universal' tube for transport back to the
- 25 laboratory.

3. DNA preparation from feather tips:

- Cut off and save the proximal 1 cm of the feather (ie the non-barbed part which is attached to the skin and which contains the pulp – see
- 30 photographs). Discard the distal barbed part of the feather.

For each chicken, place the 8 – 10 saved feather ends in a 1.5 ml snap-cap Eppendorf tube.

- Add 500 μ l proteinase K sample buffer (TNE buffer containing 0.5% SDS) containing 100 μ g proteinase K (add proteinase K just before use).
- 5 • Incubate at 50°C in a water-bath for 1.5 – 2 hours.
- Microcentrifuge the tubes (6000 rpm, 10 min), to 'pellet' feather tips and debris.
- Transfer supernatant to a new snap-cap tube (if the feathers were from brown birds, the supernatants will be brown due to the presence of melanin).
- 10 • Add an equal volume (500 μ l) of phenol to the supernatant
- Vortex
- Centrifuge at 13000 rpm, 2 min
- Transfer the upper (aqueous) phase to a new snap-cap tube
- 15 • Add an equal volume (500 μ l) of phenol to the supernatant again.
- Vortex
- Centrifuge at 13000 rpm, 2 min.
- Transfer the upper phase to a new snap-cap tube.
- Add an equal volume (500 μ l) of cold chloroform.
- 20 • Vortex.
- Centrifuge at 13000 rpm, 2 min.
- Transfer the upper phase to a 1.5 ml screw-cap tube.
- Add 1 ml filtered 100% ethanol.
- Add 50 μ l of 3M Sodium Acetate.
- 25 • Gently mix by inverting the tube, and leave at room temperature for 20 minutes (the DNA will become visible as it precipitates).

Centrifuge at 13000 rpm, 2 minutes, to pellet the DNA (if white chickens were used, pellet will be white; if brown chickens used, pellet brown).

- Discard the supernatant.
- 5 • Rinse the pellet twice with 500 μ l of 70% cold ethanol, by gently running the ethanol down the side of the tube, then pouring off (take care not to dislodge the DNA pellet).
- Cover the open top of the tube with Parafilm, and make several puncture holes in the Parafilm using a needle.
- 10 • Place tubes in a vacuum drier for about 5 minutes to dry the pellet (alternatively air-dry).
- Re-suspend the pellet in 50 μ l PCR quality water by gentle vortexing.
- Determine the concentration of the DNA preparation using a spectrophotometer.
- 15 • Adjust the concentration to 1 mg/ml in water.
- Store at -20°C

4. 132 bp repeat PCR:

- Set up reactions on ice, in 0.5 ml Eppendorf tubes, in a PCR-dedicated
20 cabinet, using PCR-dedicated pipettes and autoclaved tips.
- Prepare a 'master mix' containing the following reagents for the appropriate number of samples (volumes given per reaction):

25

30

Component	Volume per 20 μ l reaction	Final Concentration
Forward primer (10 μ M)	1.0 μ l	0.5 μ M
Reverse primer (10 μ M)	1.0 μ l	0.5 μ M
10 x Taq buffer	2.0 μ l	
MgCl ₂ (25 mM)	1.6 μ l	2 mM
DNTP mix (10 mM)	0.5 μ l	0.25 mM
Taq gold DNA polymerase	0.1 μ l	25 units/ml
Water containing 800 μ g/ml BSA*	12.8 μ l	10 μ g BSA/reaction

- Vortex to ensure complete mixing.
 - Aliquot 19 μ l 'master mix' into autoclaved 0.5 ml snap-cap Eppendorf tubes.
 - 5 • Add 1 μ g DNA (= 1 μ l of 1 mg/ml preparation, or larger volume if DNA preparation less concentrated).
 - Vortex to ensure complete mixing.
 - (Our thermal cycler has a heated lid, so we do not need to overlay the reactions with mineral oil).
 - 10 • Run on a thermal cycler using the following cycling parameters:
- | | | |
|------|----------|-------------|
| 95°C | 2 min | 1 cycle |
| 95°C | 1 min) | |
| 50°C | 30 sec) | x 40 cycles |
| 72°C | 1 min) | |
| 72°C | 10 min | 1 cycle |
- Analyse reaction products on an agarose gel.

The results of the above experiment are shown in Figure 7.

Notes

The sense primer is located 65 bp upstream of the repeat and the antisense primer is located 105 bp downstream of the repeat. The expected band size
5 is therefore 302 bp for a single repeat (i.e. $65 + 132 + 105$), 434 bp for a double repeat ($302 + 132$) and 566 bp for a triple repeat ($434 + 132$) etc. Rispens vaccine strain produces many tandem repeats.

CLAIMS

1. A method of detecting a virus in an avian tissue sample comprising:
5 extracting genetic material from an avian tissue sample; and testing the extracted genetic material to detect any genetic material from the virus; characterised in that the avian tissue sample is derived from one or more feathers of the axillary tract.
- 10 2. A method of detecting a virus as claimed in Claim 1 wherein the method provides quantitative information on the amount of the virus in the sample.
3. A method of detecting a virus as claimed in Claim 1 or 2 wherein the
15 method is specific for MDV.
4. A method as claimed in Claim 3 wherein the method is specific for MDV serotype 1.
- 20 5. A method of detecting MDV as claimed in Claim 4 wherein the method is specific for MDV-1 Rispens strain CVI 988.
6. A method as claimed in Claim 5 wherein the method comprises:
25 (i) providing forward and reverse primers for a nucleic acid polymerase, which primers are selected from the nucleotide sequence which flanks the 132 bp repeat nucleotide sequence of MDV;
(ii) amplifying nucleic acid sequences between the primers;

- (iii) detecting the number of 132 bp repeat sequences in the amplified nucleic acid sequences; and
- (iv) relating the number of 132 bp repeat sequences to the identity of the viral nucleic acid and thereby identifying the type of MDV in the tissue sample.

7. A method as claimed in any one of Claims 1 to 5 which comprises

- (a) providing a polynucleotide sequence which is capable of binding specifically to a virus-specific target polynucleotide;
- (b) contacting the extracted genetic material with a probe whereby the probe binds specifically to its target viral polynucleotide;
- (c) determining whether the probe has bound to its target viral polynucleotide; and
- (d) determining whether the sample contains the virus on the basis that the presence of the target polynucleotide indicates the presence of the virus in the sample.

8. A method as claimed in Claim 7 wherein the step (c) of determining whether the probe has bound to a target polynucleotide comprises amplifying a region of the target polynucleotide, which region comprises the binding site of the probe.

9. A method as claimed in Claim 8 wherein amplification is primed by the following primers:

Forward primer (GGT CTG GTG GTT TCC AGG TGA)

Reverse primer (GCA TAG ACG ATG TGC TGC TGA)

10. A method as claimed in Claim 9 wherein the probe has the sequence
5' AGA CCC TGA TGA TCC GCA TTG CGA CT 3'
- 5 11. A method as claimed in any one of Claims 7 to 10 wherein the probe
is labelled fluorescently and wherein the step of determining whether
the probe has bound to a target polynucleotide comprises
determining the fluorescent emissions of the probe.
- 10 12. A method of detecting a virus as claimed in any preceding claim
wherein the method involves the use of a PCR reaction
13. A method as claimed in Claim 12 wherein before said PCR reaction
is carried out, the extracted genetic material to be tested is treated
15 with an agent to overcome the inhibitory effect of any feather tissue
factor which may be present.
14. A method of detecting a virus as claimed in Claim 13 wherein the
agent is selected from one or more of bovine serum albumin; porcine
20 (pig) albumin; and ovine (sheep) albumin.
15. A method as claimed in any preceding claim further comprising a
step of furnishing the results of the method in an intelligible format.
- 25 16. A method as claimed in Claim 15 wherein the results are furnished
by recording or storing the results of the method on an information
carrier.

17. An isolated avian tissue sample from one or more feathers from the axillary tract.

5 18. An isolated avian tissue sample as claimed in Claim 17 wherein the proximal portion of the axillary tract feather is isolated from the distal (barbed) portion of the feather.

19. A genetic material — containing extract from an isolated avian tissue sample as claimed in Claim 17 or 18.

10 20. A sample or extract as claimed in Claim 17, 18 or 19 wherein the sample is stored in a form suitable for transport to a separate location.

15 21. A sample as claimed in Claim 20 wherein the sample or extract is contained in a sealed container.

ABSTRACT

Assay Methods

5

The invention relates to methods of detecting a virus in an avian tissue sample wherein genetic material derived from feathers is tested for the presence of genetic material from the virus.

10

FIGURE 1(a).

Figure 1(a)




Mature natal down

1st Generation


Early immature feather pushing out down, now postmature.



Midimmature feather starting to emerge from sheath

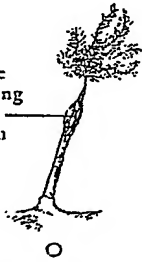

Down feather has fallen off
Late immature feather emerging from sheath
Pulp still present in ensheathed portion

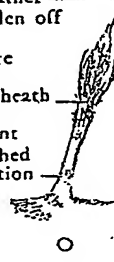

Mature 2nd generation feather emerged from sheath



Pulp entirely resorbed

2nd Generation

Early immature 3rd generation feather pushing out 2nd generation feather, now postmature



Late immature feather

Midimmature stage not shown




Mature 3rd generation feather



3rd Generation

4th Generation

FIGURE 1 (b)

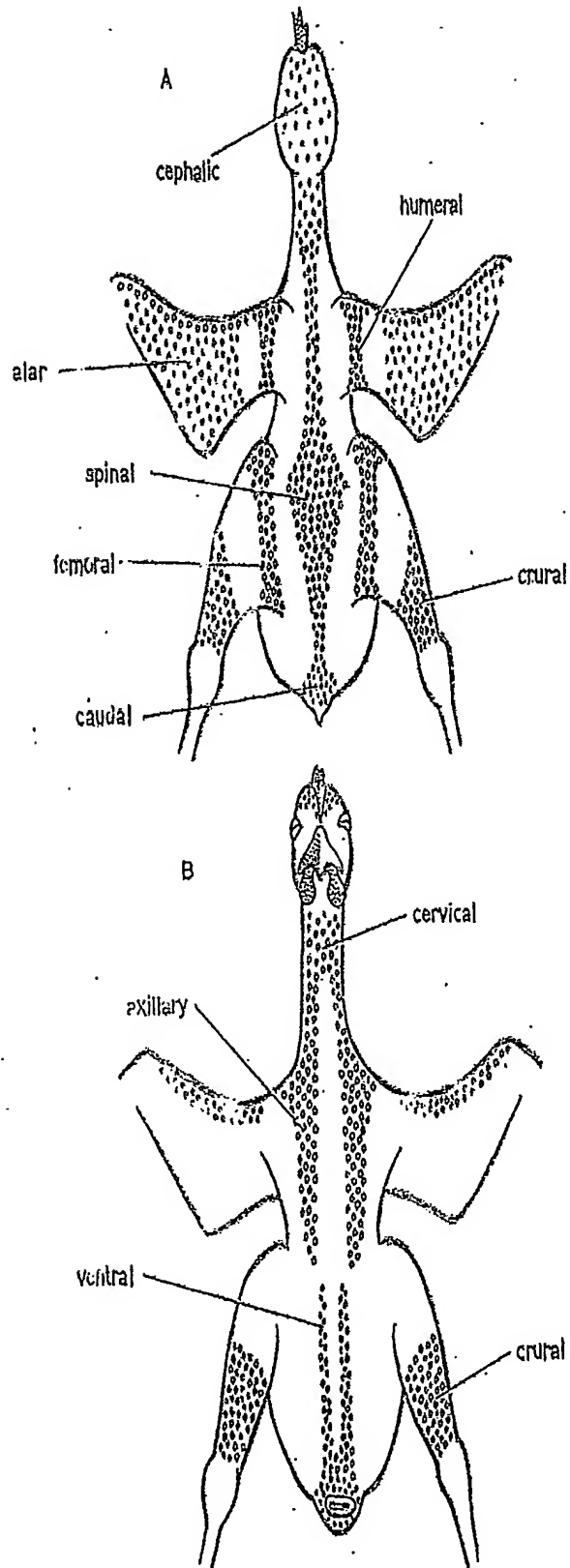

~~Fig. 1~~ The feather tracts. A. Dorsal, B. Ventral

FIGURE 2(a)

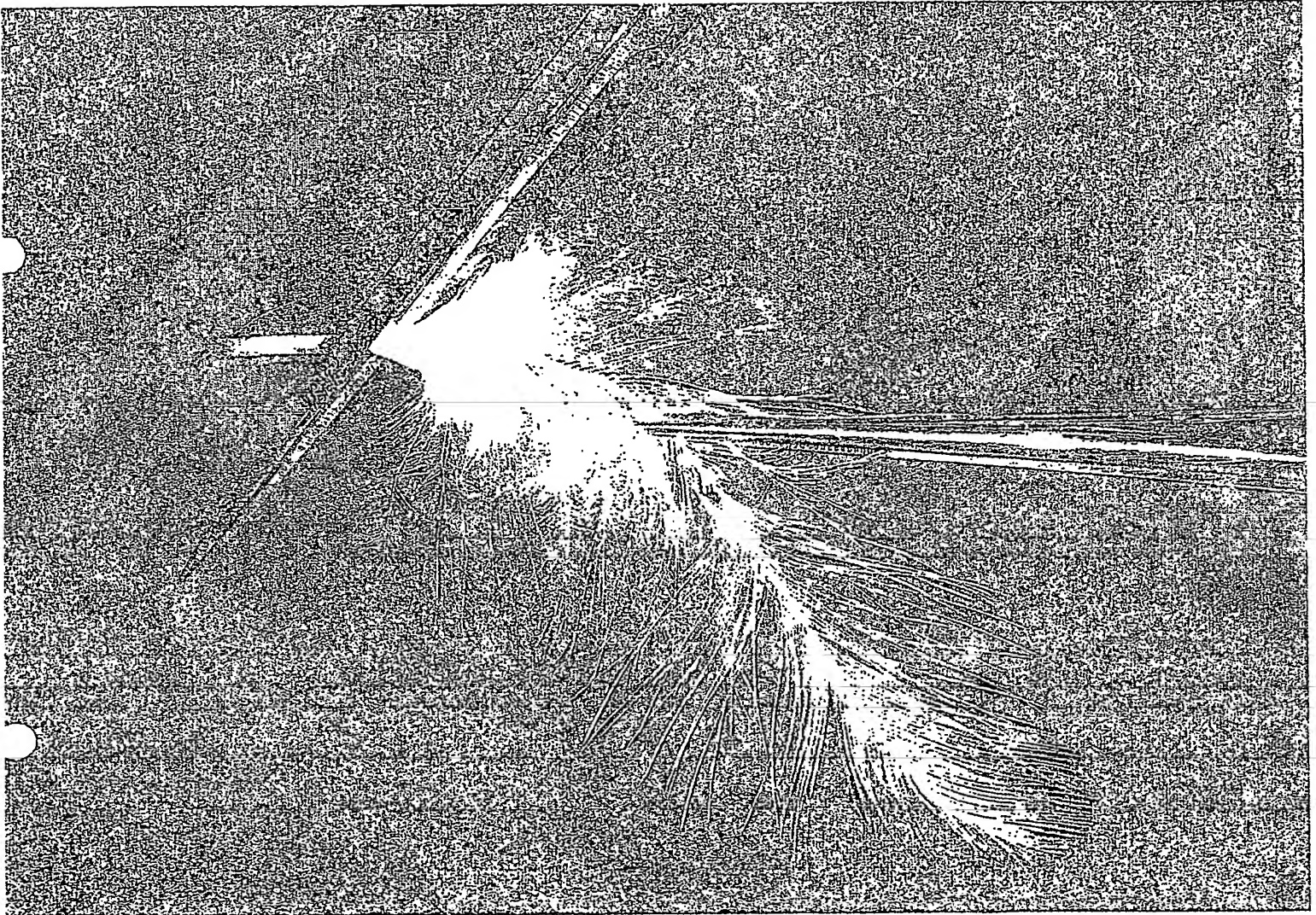


FIGURE 2(b)




5/13

FIGURE 2(c)



6/13

FIGURE 2(d)

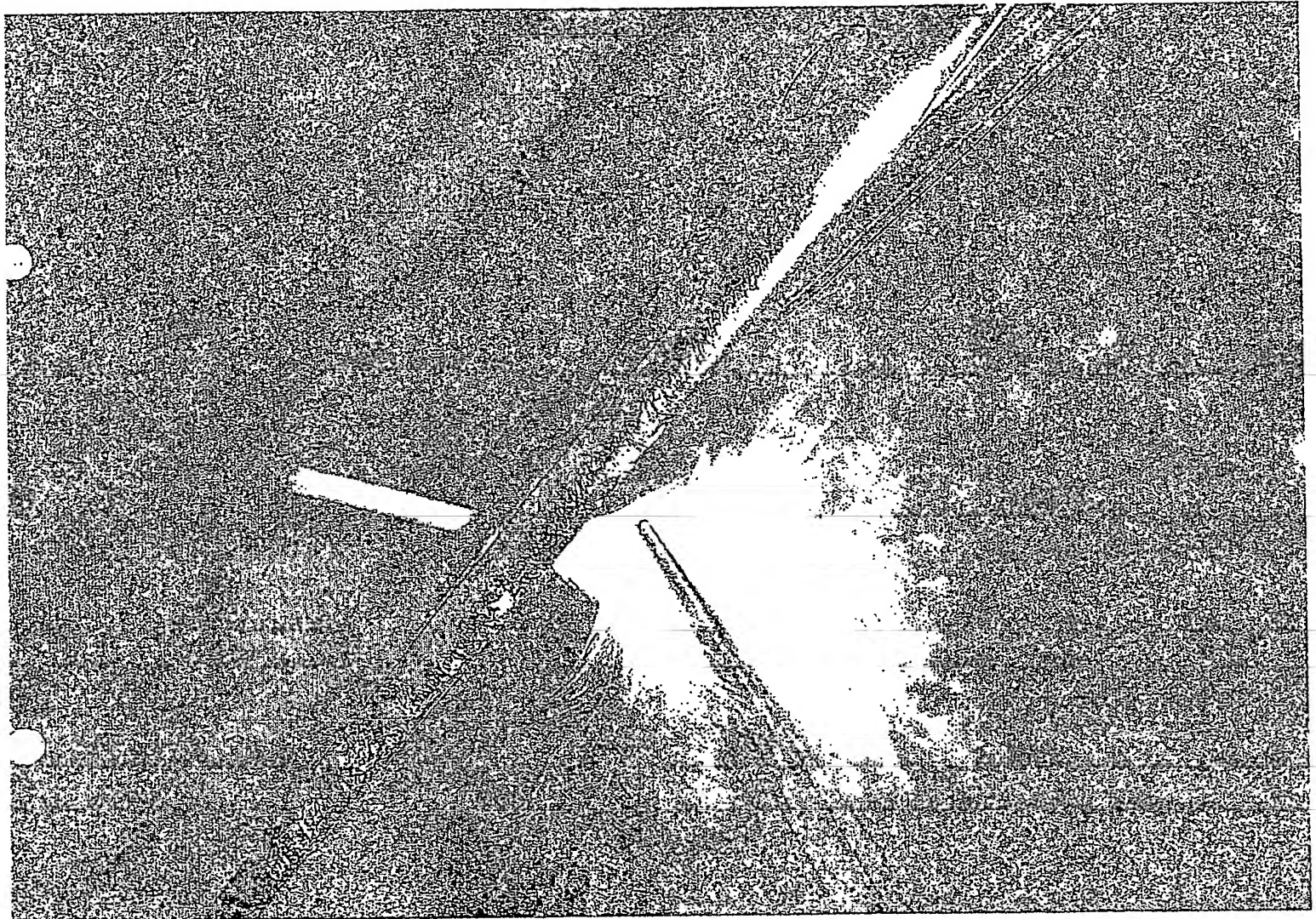

7/13

FIGURE 3

Agarose gel

Validation of technique

1. Data response of serotype 1 CVI988 vaccine as detected by ineq primers (78bp product) with Taqman technique. CVI 988 vaccine was harvested from a bacterial artificial chromosome (BAC) vector.

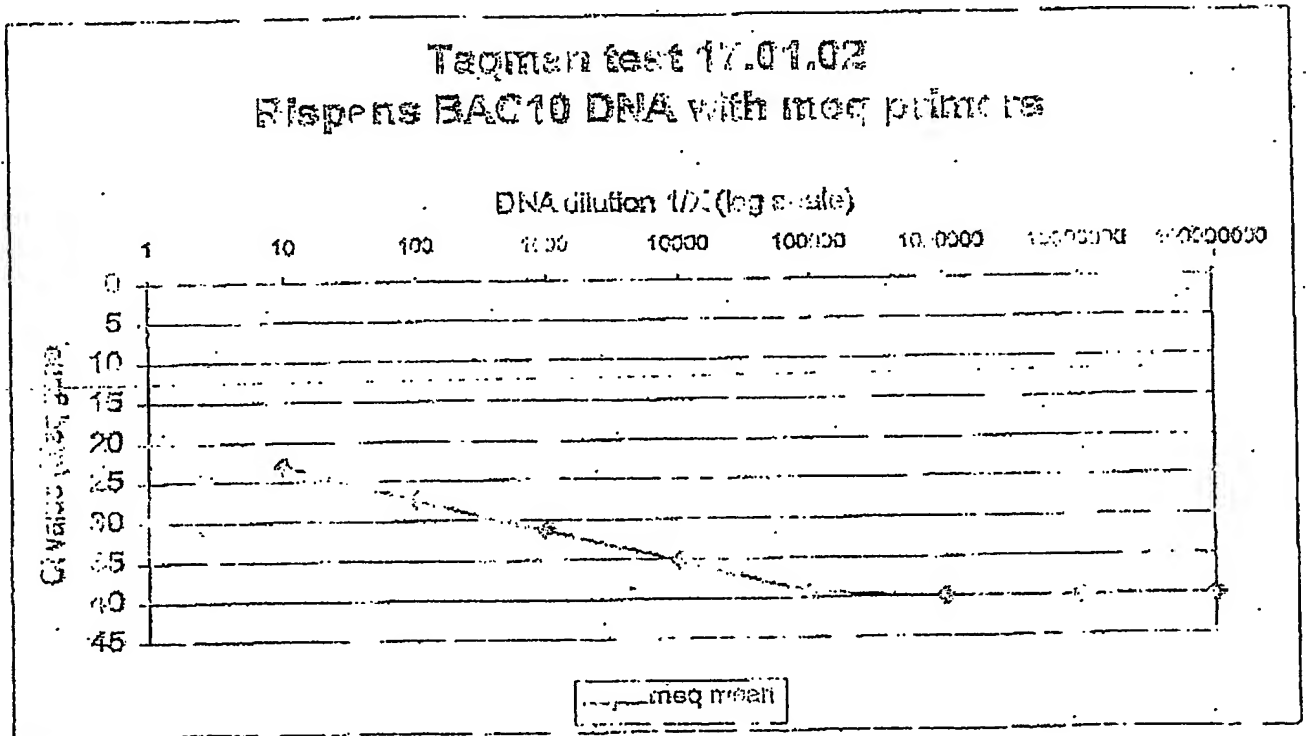


FIGURE 4

Dose response of CVI 988 vaccine in chicken embryo fibroblasts (CEF) and bacterial artificial chromosomes (BAC)


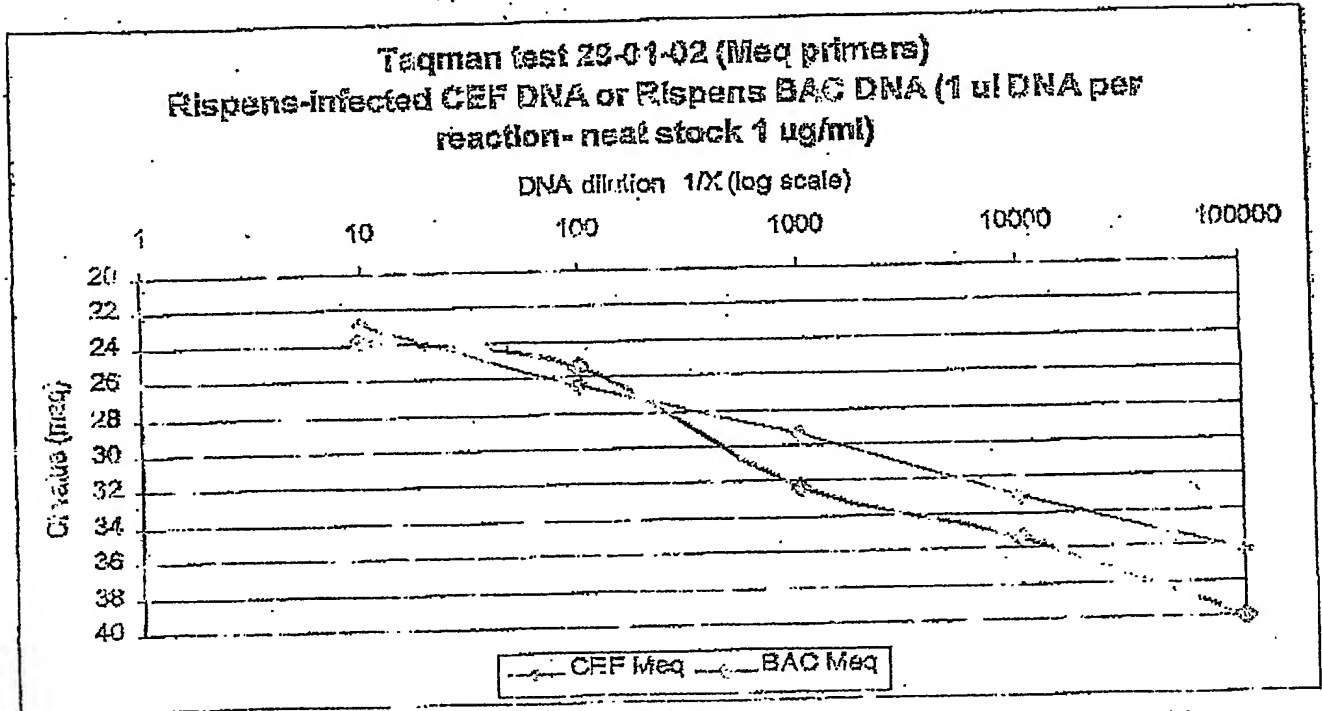


FIGURE 5

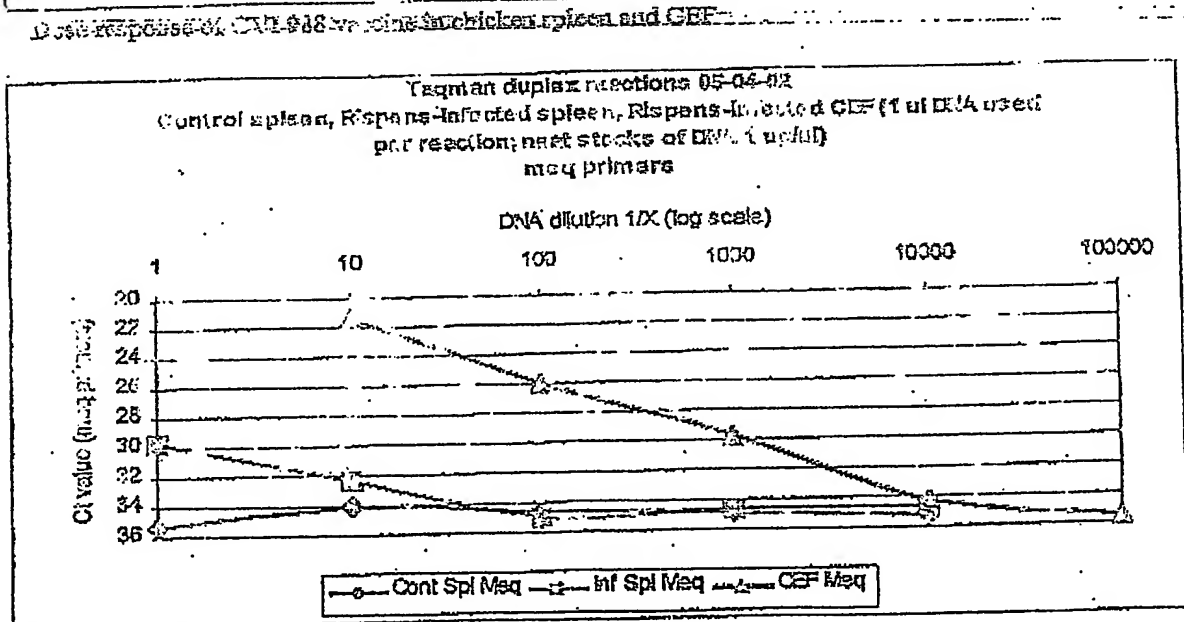

9/13

FIGURE 6

Time-course of replication of CVI988 vaccine in chicken feathers demonstrating peak replication 15-20 days post infection.

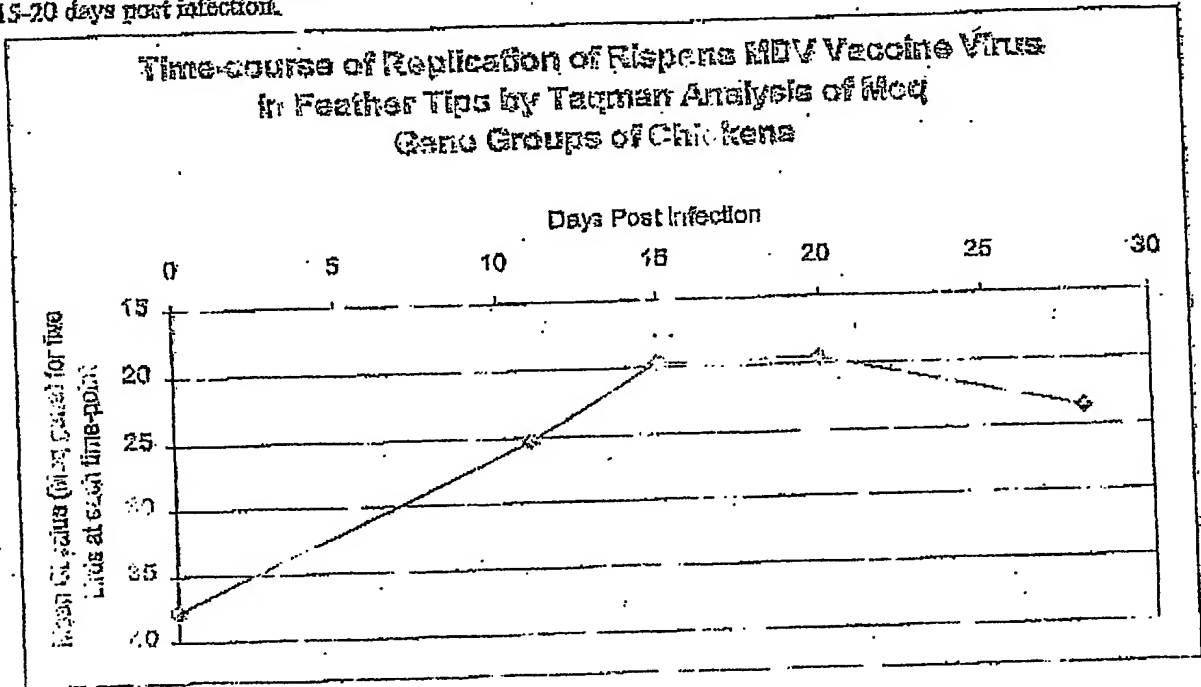


Figure 7A

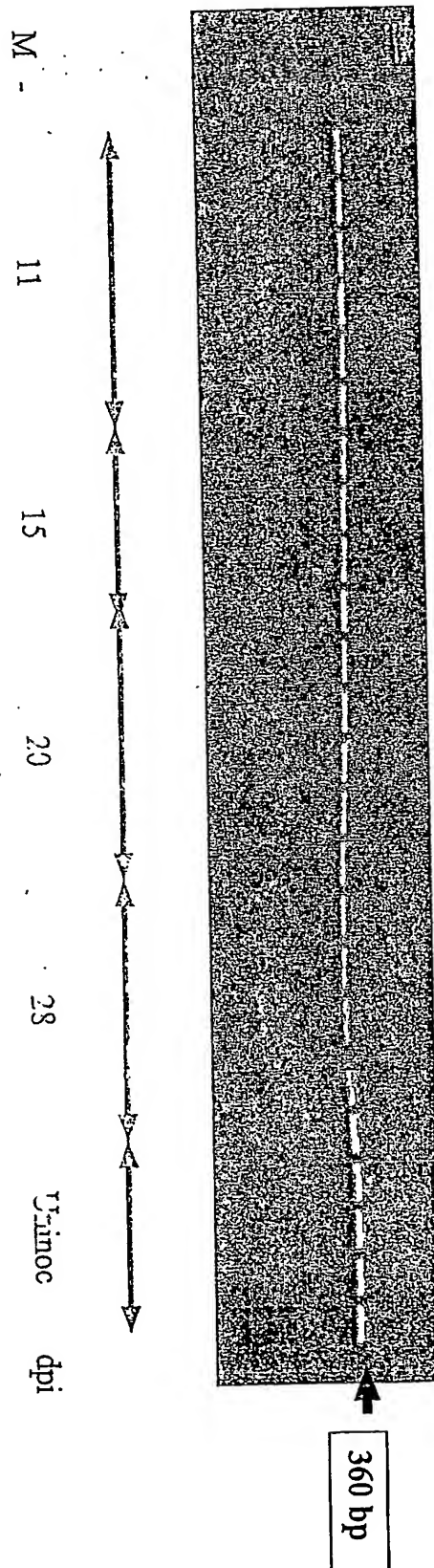


Figure 7B

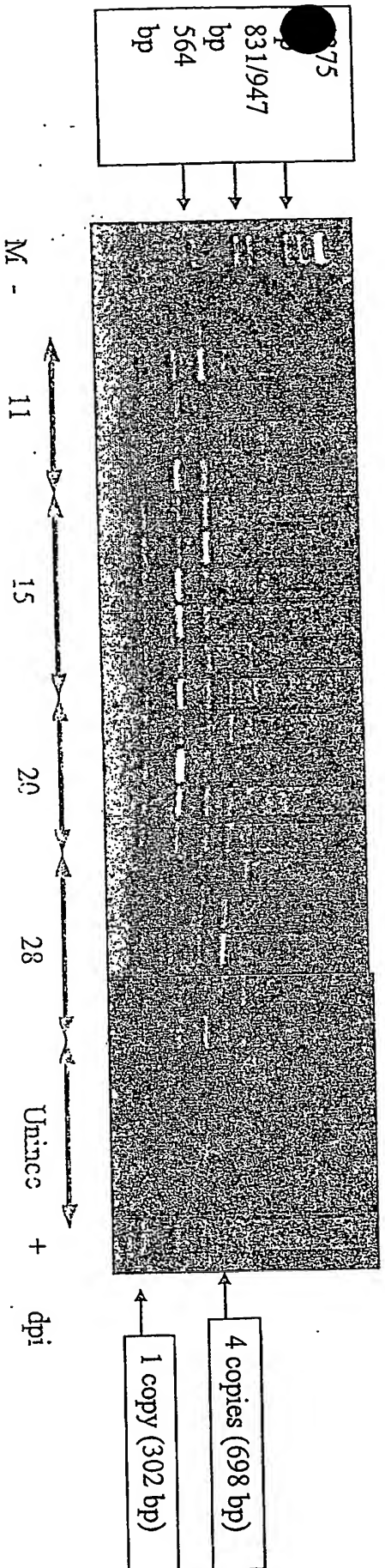


Figure 8

ID HEHSMDVE standard; DNA; VRL; 2466 BP.
 XX
 AC M89471; L09062;
 XX
 SV M89471.1
 XX
 DT 29-MAY-1992 (Rel. 32, Created)
 DT 04-MAR-2000 (Rel. 63, Last updated, Version 7)
 XX
 DE Gallid herpesvirus type 1 40kd MDV Eco Q protein gene, complete cds.
 XX
 KW leucine zipper protein; MDV Eco Q protein.
 XX
 OS Gallid herpesvirus 1
 OC Viruses; dsDNA viruses, no RNA stage; Herpesviridae; Alphaherpesvirinae;
 OC Infectious laryngotracheitis-like viruses.
 XX
 FN [1]
 RC Erratum: [Proc Natl Acad Sci U S A 1993 Mar 15;90(6):2556]
 RP 1-2466
 RX MEDLINE; 92237304.
 RA Jones D., Lee L., Liu J.L., Kung H.J., Tillotson J.K.;
 RT "Marek disease virus encodes a basic-leucine zipper gene resembling the
 RT fos/jun oncogenes that is highly expressed in lymphoblastoid tumors";
 RL Proc. Natl. Acad. Sci. U.S.A. 89(9):4042-4046(1992).
 XX
 DR SPTREMBL; Q69306; Q69306.
 XX


Key	Location/Qualifiers
source	1..2466
	/db_xref="taxon:10386"
	/organism="Gallid herpesvirus 1"
	/strain="GA"
	/tissue_lib="pBr328"
enhancer	25..38
	/note="CREB site and extended palindrome"
variation	88..89
	/note="strain differences"
	/replace="ggg"
TATA_signal	146..152
GC_signal	329..334
variation	340
	/replace="c"
CDS	355..1374
	/codon_start=1
	/db_xref="SPTREMBL:Q69306"
	/note="40 kd"
	/standard_name="~40 kD MDV Eco Q protein (meq)"
	/product="MDV Eco Q protein"
	/protein_id="AAB48631.1"
	/translation="MSQEPEPGAMPYSPADDPSPDLDSLSTSRKKRKSHDIPNSPSK HPFPDGLSEEEKQKLERRRRKRNDAARRRRRKQTDYVDKLHEACEELQRANEHLRKEIR DLRTCTSLRVQLACHEPVCMAVPLTVTLGLLTPHPDVPEPPICTPPPPSPDEPNAP HCSGSQPPICTPPPPDTEELCAQLCSTPPPPISTPHIIYAPGPSPLQPPICTPPPPDAE ELCAQLCSTPPPPICTPHSLFCPPQPPSPGIFPALCPVTEPCTPPSPGTVYAQLCPVG QAPLFTSPPHPAPEPERLYARLTEDPEQDSLYSQIYIQFPSTQSTVWWFPGDGRP"
variation	423
	/replace="c"
variation	458
	/replace="g"
variation	486

FT /replace="a"
 FT variation 583
 FT /replace="g"
 FT variation 1113
 FT /replace="c"
 FT misc_recomb 1353
 FT /organism="GA"
 FT variation 1387
 FT /replace="t"
 FT variation 1481
 FT /note="differences in GA isolates"
 FT /replace="t"
 XX
 SQ

Sequence 2466 BP; 610 A; 639 C; 586 G; 631 T; 0 other;
 gaattcgggtg atataaagac gatagtcattg catgacgtgg ggggctggat cgactgatata 60
 ctaatgggttc gggagtgata cggagacggg gggggggggg aaatgatcga tttataccta 120
 cctcttaaat aaactattgc tcctttataa aatgacaggt gaattgtgac cgttcgcgaa 180
 cgtgtaattc ttcaataactt tcgggtctgt ggggtgtgct tttttaatta ttattttggt 240
 tcggggagggt tgggtgctgga atgttaagaa taaattccgc acactgattc ctaggcaggc 300
 gtctcttgca ggtgtatacc agggagaagg cgggcacggg acagggtgaa agagatgtct 360
 caggagacag agccgggcgc tatgccctac agtcccgtc acgatccgtc cccctcgat 420
 ctttctctcg ggtcgacttc gagacggaaa aaaaggaaaa gtcacgacat cccaacagc 480
 ccctocaaac accccttccc tgacgyccta tctgaggagg agaaacagaa gctggaaagg 540
 aggagaaaaa ggaatcgtga cgcgcctcgg agaagacgca ggaagcagac ggactatgta 600
 gacaaactcc atgaagcatg tgaagagctg cagagggcca atgaacacct acgtaaggaa 660
 attcgagatc taaggactga gtgcacgtcc ctgcgtgtac agttggcttg tcatgagcca 720
 gtttgcccta tggcgggtacc cctaacgggtg acccttggac tgcttaccac cccgcacgat 780
 cccgttcctg aacctcccat ttgcactcct ccacctccct caccggatga acctaacgt 840
 ccacattgct ccggttccca acctcctatc tgtacccccc ctccctccga tacggaggaa 900
 ctttgccccc agctctgctc gacccaccca cctcccatct ctaactccca tattatctac 960
 gctccggggc cttccccctt ccaacctcct atctgtaccc cccctcctcc cgatgcggag 1020
 gagctttgctg ccagctctg ctgcacccca ccacctccca ctgtactcc ccattccctc 1080
 ttctgccctc ccagcctcc atctccggag ggaatcttc ctgcattgtg tctgtttacc 1140
 gagcgtgta cccctccatc gccggggagc gtttacgtc agctttgtcc tgttggccag 1200
 gctccctttt ttaccccatc tccccacat ccggtcccg agccggagag gctttatgct 1260
 cgtcttaccg aggatccga acaggattcc ttgtattcgg gccagattta tattcagttt 1320
 ccctcgata ctacgtctac ggtctggtgg ttccagggtg acgggagacc ctgatgatcc 1380
 gcattgcgac tctcagcagc acatcgtcta tgcccatgt ttcttctccc ctagtatat 1440
 ataatagttt tcatagtltc gggaagatca acataaagga aagggttaaa ggcattat 1500
 atcgatttac tgacataaaa aaactcctctg gggtacaaa ttttccctta ccgtgtagct 1560
 tagactcgga agaactat 1620
 cgaagtatga gataaactta gctatgtgga aaacttctgg ggcaacatct ctcgcccca 1680
 gactgcttaa atggcaaatt ctggttctat acagaacggt tggggaagg gggggggggg 1740
 glatatggag tattattcgg gatattgctt ctatgaagct gcggttaagt ttccaggctc 1800
 aaaaactatg cctggctgtt ttttttttta ctagggatat ggacatcgca cattaaggaa 1860
 tattaagat aacaggatgg acattcggat gtaaaaggaa taagcgaac ctttcgcaga 1920
 tgtgagttaa tgcagtctcg tataattcgg tgggtctgat taggttatcg taaggaaaca 1980
 cagattgat ctctcatccg cgtcccagca atcaggccta tgtccctctc ctgtggccag 2040
 ctactggct gtgcactgtg cgattctaag tgetacagtc gtgagcagat caatggatcg 2100
 gggctcgcgc acactactg taattaaata ttcgtttatg aattatgcaa atatgcacag 2160
 ataatatata cagggatgca cagacatact cctatgcacc gatacacagg cacataggca 2220
 gatgtcgaca ttaacgaata tacaggcaag gacctccagg aacatatgga aaataacctca 2280
 tcgcagagac gcttatgcag gagtaactct cgttaagtcg ttactggatt gtaacggcta 2340
 tccggagact ctcttccct tttgcttgtt cactgtgcgg cattattaca tttacaccgg 2400
 taatgctgcg catgaaagag cgaacggaac gaggtctgta cgacattaca agaatagttt 2460
 gaattc 2466

PCT Application

EP0311360

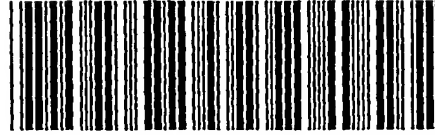

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☒ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☒ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☒ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.